



Aesculap Orthopaedics OrthoPilot®



KneeSuite™ – ACL Module

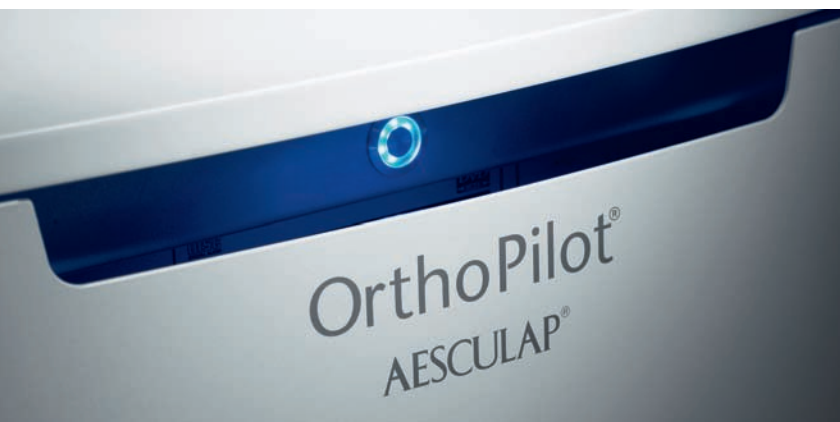

KneeSuite™ ACL Module Anterior Cruciate Ligament

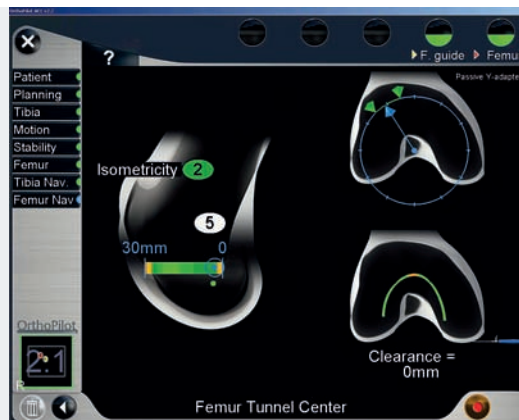

OrthoPilot
AESCULA

Powered by the OrthoPilot® Navigation System, the ACL reconstruction software module is easy to use and offers the same innovative technology and capabilities as the KneeSuite™ and HipSuite™ applications.

Precision



By referencing a few anatomic landmarks, the surgeon has the ability to determine the optimal location of the tibial and femoral bone tunnels based on the anatomic and/or isometric criteria.

Additionally, OrthoPilot® provides information about potential impingement to ensure optimal function of the graft.

Integrated Documentation


The OrthoPilot ACL module offers a pre-op and post-op stability test. The stability test is documented without exposing the patient to additional radiation. This information can be recorded, printed and placed in the patient's file.

lot[®]
P[®]

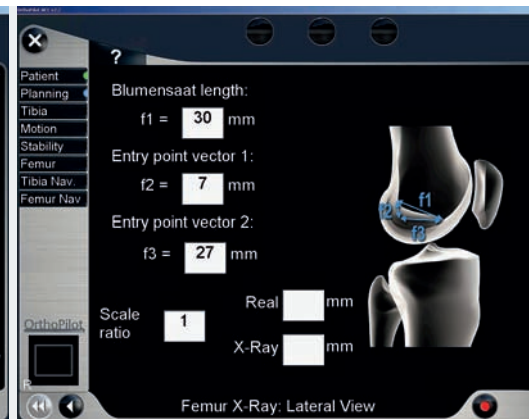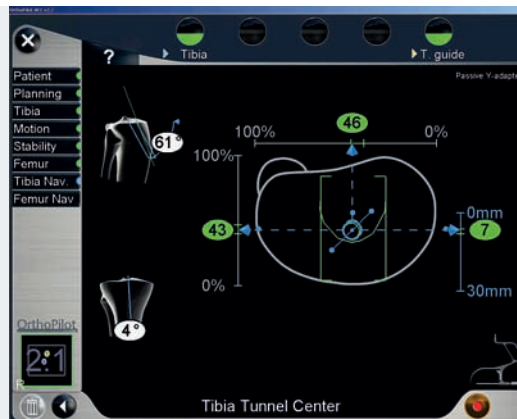

Features

- > Optimized and streamlined instruments
- > Optional pre-op planning of the femoral tunnel position can be obtained using the quadrant method on existing pre-op X-ray images
- > Reconstruction according to anatomic and/or isometric criteria
- > Display of potential impingement
- > State-of-the-art passive transmitter technology for simple and safe intraoperative handling
- > Pre-op and post-op stability tests (Lachman test & internal/external rotation)

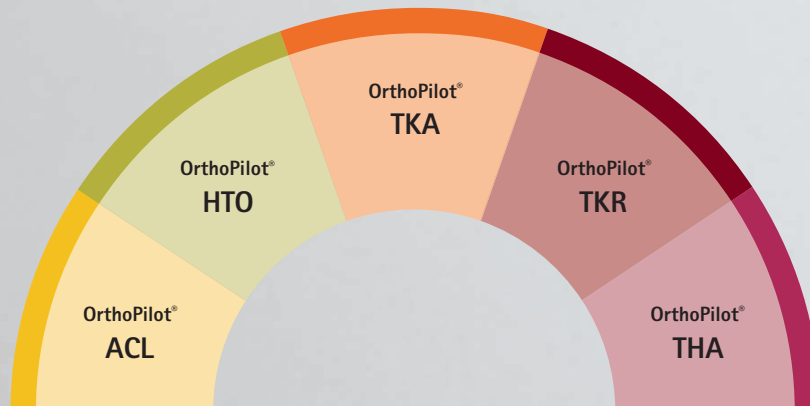

AESCULAP®

Implant Systems

Aesculap Implant Systems, Inc.
3773 Corporate Parkway
Center Valley, PA 18034

Phone: 866-229-3002
www.aesculapimplantsystems.com

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogs and price lists and to take legal actions.

©2008 AESCULAP. ALL RIGHTS RESERVED. PRINTED IN THE USA. DOC579 REV. B 1M 10/08
Aesculap is an equal opportunity employer